

ISSN:2348-2079

Volume-8 Issue-3

International Journal of Intellectual Advancements and Research in Engineering Computations

A validated RP-HPLC method for the simultaneous estimation of cephlexin and kanamycin in its pure and pharmaceutical dosage form

Dr Gutta Subbarao and M.Supriya

St. Xavier Institute of Pharmacy, Deenapur, Phirangipuram, Guntur, India.

ABSTRACT

A new method was established for simultaneous estimation of Kanamycin and Cephalexin by RP-HPLC method. The chromatographic conditions were successfully developed for the separation of Kanamycin and Cephalexin by using symmetry C18 column $(4.6\times150\text{mm})$ 5 μ , flow rate was 1ml/min, mobile phase ratio was (70:30 v/v) methanol: phosphate buffer $(KH_2PO_4 \text{ and } K_2HPO_4)$ pH 3 (pH was adjusted with orthophosphoric acid), detection wave length was 258nm. The instrument used was WATERS HPLC Auto Sampler, Separation module 2695, photo diode array detector 996, Empower-software version-2. The retention times were found to be 2.403 mins and 3.907 mins. The % purity of Kanamycin and Cephalexin was found to be 100.27% and 99.87% respectively. The system suitability parameters for Kanamycin and Cephalexin such as theoretical plates and tailing factor were found to be 2294, 1.27 and 4891 and 1.03, the resolution was found to be 8.67. The analytical method was validated according to ICH guidelines (ICH, Q2 (R1)). The linearity study n Kanamycin and Cephalexin was found in concentration range of $50\mu g$ -250 μ g and $5\mu g$ -50 μ g and correlation coefficient (μ) was found to be 0.999 and 0.999, recovery was found to be 99.56% and 99.48%, %RSD for repeatability was 0.27 and 0.40, % RSD for intermediate precision was 0.27 and 0.94 respectively. The precision study was precise, robust, and repeatable.LOD value was 2.17 and 6.60, and LOQ value was 0.032 and 0.1125 respectively.

Keywords: RP-HPLC Method, Cephlexin and Kanamycin

INTRODUCTION

A drug is a substance which may have medicinal, intoxicating, performance enhancing or other effects when taken or put into a human body or the body of another animal and is not considered a food or exclusively a food. What is considered a drug rather than a food varies between cultures, and distinctions between drugs and foods and between kinds of drug are enshrined in laws which vary between jurisdictions and aim to restrict or prevent drug use. Even within a jurisdiction, however, the status of a substance may be uncertain or contested with respect to both whether it is a drug and how it should be classified if at all. There is no single, precise definition, as there are different meanings in drug control law, government regulations, medicine, and Colloquial usage.

In pharmacology, a drug is "a chemical substance used in the treatment, cure, prevention, or diagnosis of disease or used to otherwise enhance physical or mental well-being. Drugs may be prescribed for a limited duration, or on a regular basis for chronic disorders.

The molecules of drugs are complex, and most of them consist of many hydrogen and carbon atoms, a few oxygen atoms, and one or a few nitrogen atoms. Drugs may also have no nitrogen atoms in it and many may have chlorine atoms in it, such as chloral hydrate.

Recreational drugs are chemical substances that affect the central nervous system, such as opioids or hallucinogens. They may be used for perceived beneficial effects on perception, consciousness, personality, and behaviour. Some drugs can cause addiction and/or habituation.

Analytical chemistry is applied throughout industry, medicine, and all the sciences. Consider a few examples. The concentration of oxygen and of oxygen and of carbon dioxide are determined in millions of blood samples every day and used to diagnose and treat illnesses. Quantities of hydrocarbons, nitrogen oxides, and carbon monoxide present in automobile exhaust gases are measured to assess the effectiveness of smog-control devices. Quantitative measurements of ionized calcium in blood serum help diagnose parathyroid disease in humans. Quantitative determination of nitrogen in foods establishes their protein content and thus their nutritional value. Analysis of steel during its production permits adjustment in the concentration of such elements as carbon, nickel, and ductility. The mercaptan content of household gas supplies is monitored continually to ensure that the gas has a sufficiently noxious odour to warm of dangerous leaks. Farmers tailor fertilization and irrigation schedules to meet changing plant needs during the growing season. Gauging these needs from quantitative analyses of the plants and the soil in which they grow.

Author for correspondence:

Quantitative analytical measurements also play a vital role in many research areas in chemistry, biochemistry, biology, geology, physics, and the other sciences.

For example, quantitative analytical measurement of potassium, calcium, and sodium ions in the body fluids of animals permit physiologists to study the role of these ions in nerve signal conduction as well as muscle contraction and relaxation. Chemists unravel the mechanisms of chemical reactions through reaction rate studies, the rate of consumption of reactants or formation of products in a chemical reaction can be calculated from quantitative measurements made at equal time intervals. Material scientists rely heavily on quantitative analysis of crystalline germanium and silicon in their studies of semiconductor device.

Quantitative analytical methods

We compute the results of a typical quantitative analysis from two measurements. One is the mass or the volume of sample being analyzed. The second is the measurement of some quantity that is proportional to the amount of analyte in the sample, such as mass, volume, intensity of light, or electrical charge. This second measurement usually completes the analysis, and we classify analytical methods according to the nature of this final measurement. Gravimetric methods determine the mass of the analyte or some compound chemically related to it. In a volumetric method, the volume of a solution containing sufficient reagent to react completely with the analyte is measured. Electro analytical methods involve the measurement of such electrical properties as potential, current, resistance, and quantity of electrical charge. Spectroscopic methods are based on measurement of the interaction between electromagnetic radiations by analytes. Finally, a group of miscellaneous methods includes the measurement of such quantities as mass-to - charge ratio of molecules by mass spectrometry, rate of radioactive decay, heat of reaction, rate of reaction, sample thermal conductivity, optical activity, and refractive index.

Pharmaceutical Analysis plays a very vital role in the quality assurance and quality control of bulk drugs and their formulations. (2) Pharmaceutical analysis is a specialized branch of analytical chemistry, which involves separating, identifying and determining the relative amounts of components in a sample of matter. Pharmaceutical analysis derives its principles from various branches of sciences like physics, microbiology, nuclear science and electronics etc. Qualitative analysis reveals the chemical identity of the sample. Quantitative analysis establishes the relative amount of one or more of these species or analytes in numerical terms nearly, any physical property or characteristic of a particular element or compound can be made the basis of a method for its analytical determination.

E.g.: Spectroscopic technique involves the absorption/ emission of radiant energy in all regions of the electromagnetic spectrum.

In recent years, several analytical techniques have been evolved that combine two or more methods into one called "hyphenated" technique e.g.: GC/MS, LC/MS etc.

The complete analysis of a substance consists of four main steps

1. Sample preparation / Sampling.

- Dissolution of the sample, conversion of the analyte into a form suitable for measurement.
- Measurement.
- 4. Calculation and interpretation of the measurement.

MATERIALS AND METHODS

Analytical method validation

Validation parameters

- Specificity
- Linearity
- Range
- Accuracy
- Precision

1. Specificity

The system suitability for specificity was carried out to determine whether there is any interference of any impurities in retention time of analytical peak. The specificity was performed by injecting blank.

2. Linearity

Preparation of stock solution

10 mg of kanamycin and 1mg of cephalexin working standard were accurately weighed and were transferred into a 10ml clean dry volumetric flask, add about 2ml of diluent and sonicate to dissolve it completely and make volume up to the mark with the same solvent.

Preparation of Level – I (50ppm of kanamycin and 5 ppm of cephalexin)

0.5 ml of stock solution was taken in to 10ml of volumetric flask and diluted up to the mark with diluent.

Preparation of Level – II (100ppm of kanamycin and 10ppm of cephalexin)

1 ml of stock solution was taken in to 10ml of volumetric flask and diluted up to the mark with diluent.

Preparation of Level - III (150ppm of kanamycin and 15ppm of cephalexin)

1.5 ml of stock solution was taken in to 10ml of volumetric flask and diluted up to the mark with diluent.

Preparation of Level – IV (200 ppm of kanamycin and 20ppm of cephalexin)

2 ml of stock solution was taken in to 10ml of volumetric flask and diluted up to the mark with diluent.

Preparation of Level -V (250 ppm of kanamycin and 25ppm of cephalexin)

2.5 ml of stock solution was taken in to 10ml of volumetric flask and diluted up to the mark with diluent.

Procedure

Each level was injected into the chromatographic system and peak area was measured. Plot a graph of peak area

versus concentration (on X-axis concentration and on Y-axis Peak area) and the correlation coefficient was calculated. Acceptance criteria.

Correlation coefficient should be not less than 0.999.

3. Range

Based on precision, linearity and accuracy data it can be concluded that the assay method is precise, linear and accurate in the range of $50\mu g/ml\text{-}250\mu g/ml$ and $5\mu g/ml\text{-}25\mu g/ml$ of kanamycin and cephalexin respectively.

4. Accuracy

Preparation of standard stock solution

10mg of kanamycin and 1mg of cephalexin working standard were accurately weighed and transferred into a 10ml clean dry volumetric flask add about 2ml of diluent and sonicate to dissolve it completely and make volume up to the mark with the same solvent (Stock solution). Further pipette out 1 ml of the above stock solution into a 10 ml volumetric flask and was diluted up to the mark with diluent.

Preparation of sample solutions

For preparation of 50% solution (with respect to target assay concentration)

5mg of kanamycin and 0.5 mg of cephalexin working standard were accurately weighed and transferred into a 10 ml clean dry volumetric flask add about 2ml of diluent and sonicate to dissolve it completely and make volume up to the mark with the same solvent (Stock Solution). Further pipette out 10 ml of the above stock solution into a 100ml volumetric flask and was diluted up to the mark with diluent.

For preparation of 100% solution (with respect to target assay concentration)

10 mg of kanamycin and 1 mg of cephalexin working standards were accurately weighed and transferred into a 10ml clean dry volumetric flask add about 2 ml of diluent and sonicate to dissolve it completely and make volume up to the mark with the same solvent (Stock solution). Further pipette out 1ml of above stock solution into a 10 ml volumetric flask and was diluted up to the mark with diluent.

For preparation of 150% solution (with respect to target assay concentration)

15 mg of kanamycin and 2 mg of cephalexin working standards into a 10ml clean dry volumetric flask add about 2ml of diluent and sonicate to dissolve it completely and make volume up to the mark with the same solvent. Further pipette out 1ml of the above stock solution into a 10ml volumetric flask and was diluted up to the mark with diluent.

Procedure

The standard solutions of accuracy 50%, 100% and 150% were injected into chromatographic system. Calculate the amount found and amount added for kanamycin and cephalexin and calculate the individual % recovery and mean % recovery values.

Acceptance criteria

The % recovery for each level should be between 98.0 to 102.0%

RESULTS AND DISCUSSION

Validation report

Specificity

The system suitability for specificity was carried out to determine whether there is any interference of any impurities in retention time of analytical peak. The study was performed by injecting blank. The chromatograms are shown in Fig.No.34-35.

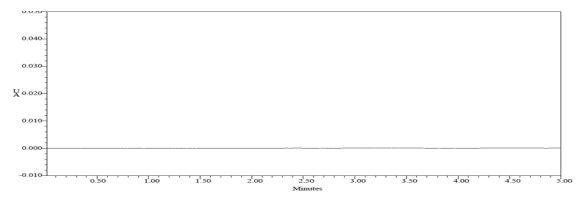


Fig 1: Chromatogram showing blank (mobile phase preparation)

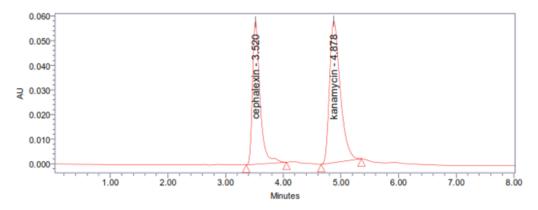


Fig 2: Chromatogram showing standard injection

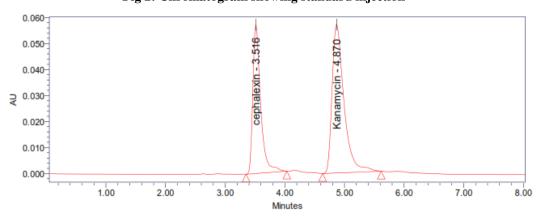



Fig 3: Chromatogram showing sample injection

The specificity test was performed for Kanamycin and cephalexin. It was found that there was no interference of impurities in retention time of analytical peak.

Linearity

The linearity study was performed for the concentration

of 50 ppm to 250 ppm and 5ppm to 25 ppm level. Each level was injected into chromatographic system. The area of each level was used for calculation of correlation coefficient. The chromatograms are shown in Fig. No.36-40 and results are tabulated in Table. No.14-15. Calibration graph for CEF and KAN are shown in Fig.No.41-42.

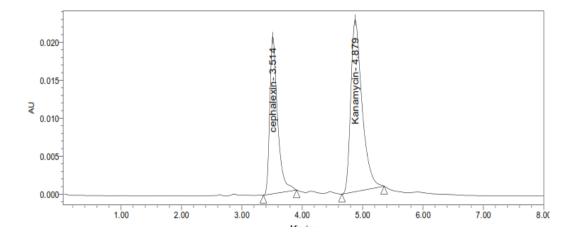


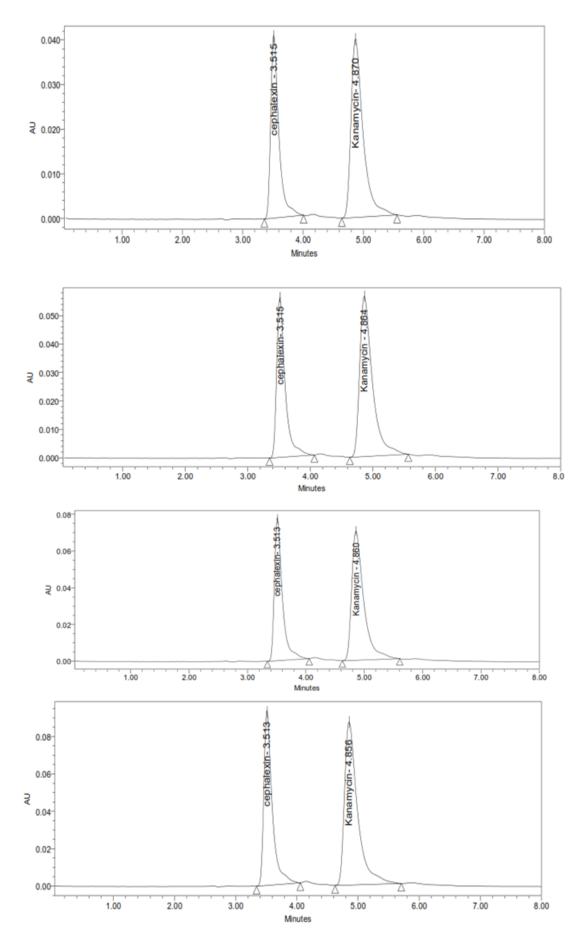


Fig 4: Chromatograms showing linearity level-1 to level 5 (50ppm-250 ppm

Table 1: Linearity Results for Cephalaxin


S.No	Linearity Level	Concentration	Area			
1	I	50 ppm	471543			
2	II	100 ppm	656277			
3	III	150 ppm	794999			
4	IV	200 ppm	946124			
5	V	250 ppm	1002139			
	Correlation Coefficient					

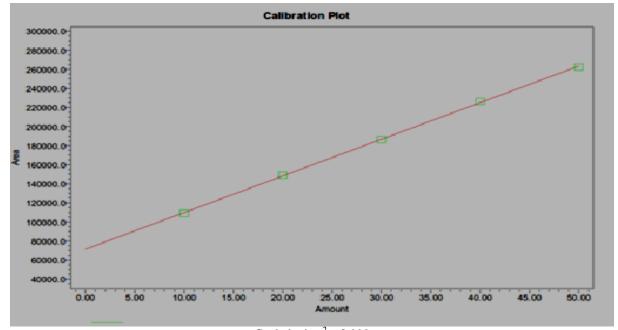

Cephalexin $r^2 = 0.999$


Table 2: Linearity Results for Kanamycin

S.No	Linearity Level	Concentration	Area
1	I	5ppm	56472
2	II	10 ppm	73841
3	III	15ppm	92655
4	IV	20ppm	111541
5	V	25ppm	130567
	0.999		

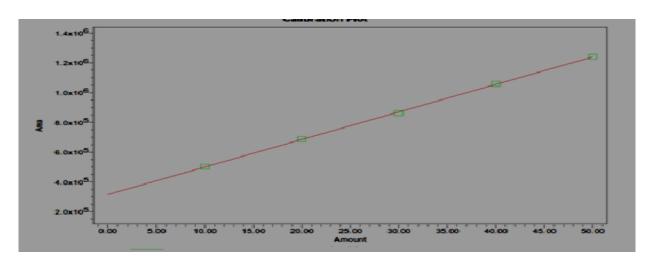

Overlay diagram

Fig 5: Showing calibration graph for Kanamycin

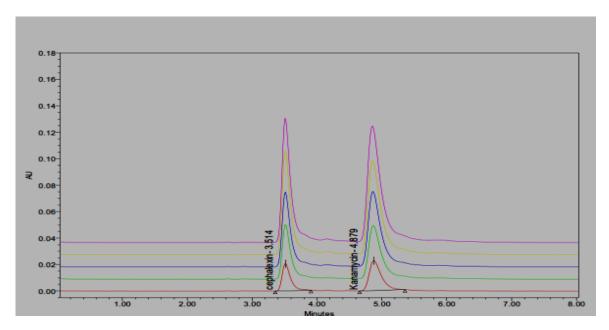


The linearity study was performed for concentration range of $50.\mu g$ - $250\mu g$ and $5\mu g$ - $50\mu g$ of kanamycin and cephalexin and the correlation coefficient was found to be 0.999 and 0.999.(NLT 0.999).

120 % for kanamycin and cephalexin. Each level was injected in triplicate into chromatographic system. The area of each level was used for calculation of % recover y. Chromatograms are shown in Fig.No.43-51 and results are tabulated in Table.No.16-17

Accuracy

The accuracy study was performed for 80%, 100% and

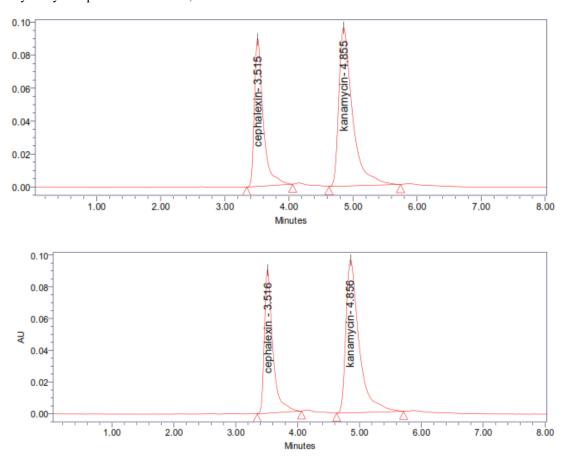


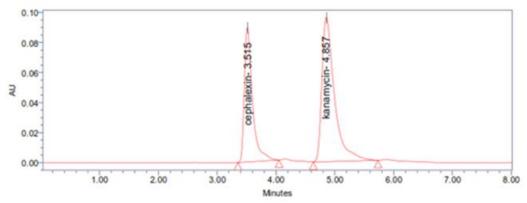



Fig 6: Chromatograms showing accuracy-80% injection-1, 2, 3

Peak Results

Tour to diffe							
Name	RT	Area	Height	USP Resolution	USP Tailing	USP Plate Count	Injection
cephalexin	3.515	852858	90288		1.79	3645	1
kanamycin	4.855	1425312	96264	4.32	1.99	2850	1
cephalexin	3.516	854502	90471		1.74	3685	2
kanamycin	4.856	1425604	96724	4.34	1.98	2844	2
Febuxostat	3.515	853960	89576		1.78	3629	3
kanamycin	4.857	1424941	96279	4.32	1.98	2847	3
	cephalexin kanamycin cephalexin kanamycin Febuxostat	cephalexin 3.515 kanamycin 4.855 cephalexin 3.516 kanamycin 4.856 Febuxostat 3.515	cephalexin 3.515 852858 kanamycin 4.855 1425312 cephalexin 3.516 854502 kanamycin 4.856 1425604 Febuxostat 3.515 853960	cephalexin 3.515 852858 90288 kanamycin 4.855 1425312 96264 cephalexin 3.516 854502 90471 kanamycin 4.856 1425604 96724 Febuxostat 3.515 853960 89576	cephalexin 3.515 852858 90288 kanamycin 4.855 1425312 96264 4.32 cephalexin 3.516 854502 90471 kanamycin 4.856 1425604 96724 4.34 Febuxostat 3.515 853960 89576	cephalexin 3.515 852858 90288 1.79 kanamycin 4.855 1425312 96264 4.32 1.99 cephalexin 3.516 854502 90471 1.74 kanamycin 4.856 1425604 96724 4.34 1.98 Febuxostat 3.515 853960 89576 1.78	cephalexin 3.515 852858 90288 1.79 3645 kanamycin 4.855 1425312 96264 4.32 1.99 2850 cephalexin 3.516 854502 90471 1.74 3685 kanamycin 4.856 1425604 96724 4.34 1.98 2844 Febuxostat 3.515 853960 89576 1.78 3629

 $\label{eq:Accuracy-100\%} Accuracy -100\% \\ \mbox{Fig 7: Chromatogram showing accuracy -100\% injection-1, 2, 3.}$

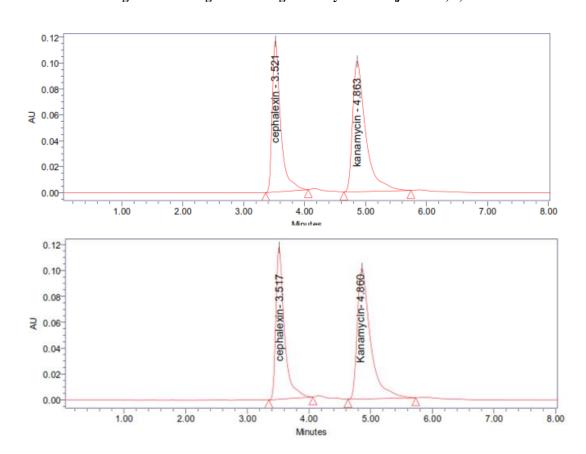

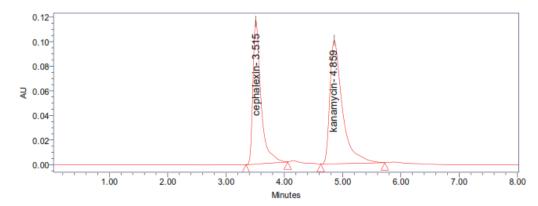

Peak Results

_								
	Name	RT	Area	Height	USP Resolution	USP Tailing	USP Plate Count	Injection
1	cephalexin	3.521	1119197	116736		1.84	3637	1
2	kanamycin	4.863	1499388	100942	4.32	2.04	2860	1
3	cephalexin	3.517	1119910	117508		1.80	3643	2
4	kanamycin	4.860	1500972	100747	4.33	2.02	2833	2
5	cephalexin	3.515	1118239	117424		1.86	3627	3
6	kanamycin	4.859	1499296	100606	4.31	2.03	2852	3

Accuracy 120%

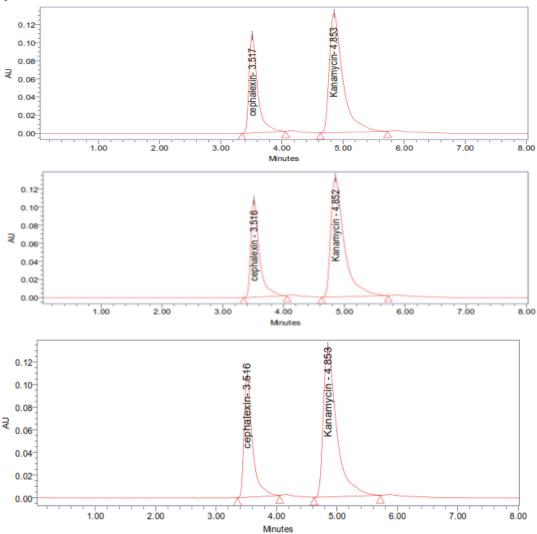


Fig 8: Chromatogram showing accuracy -120 % injection-1,2,3.

 $\label{lem:computations} \begin{tabular}{l} Copyrights @ International Journal of Intellectual Advancements and Research in Engineering Computations, \\ www.ijiarec.com \end{tabular}$

Peak Results

_								
	Name	RT	Area	Height	USP Resolution	USP Tailing	USP Plate Count	Injection
1	cephalexin	3.517	1038871	108527		1.87	3623	1
2	Kanamycin	4.853	2024658	132459	4.24	2.09	2678	1
3	cephalexin	3.516	1038553	107528		1.86	3611	2
4	Kanamycin	4.852	2021564	131567	4.22	2.09	2705	2
5	cephalexin	3.516	1037242	108249		1.89	3631	3
6	Kanamycin	4.853	2021860	131786	4.23	2.09	2694	3

Table 3: Showing accuracy results for Cephalexin

%Concentration (at specification level)	Average area	Amount added (mg)	Amount found (mg)	% Recovery	Mean recovery
50%	7371253	5	4.96	99.91%	_
100%	14634226.7	10	9.98	99.18%	99.56%
150%	2243270.7	15	15.02	99.60%	•

Table.No.4. Showing accuracy results for Kanamycin

%Concentration (at specification level)	Average Area	Amount added (mg)	Amount found (mg)	% Recovery	Mean recovery
50%	484733	0.5	0.99	99.53%	
100%	967998	1.0	1.05	99.38%	99.47%
150%	145437	1.5	1.495	99.52%	

The accuracy study was performed for % recovery of Kanamycin and cephalexin. The % recovery was found to be 99.18% and 99.91% respectively (NLT 98% and NMT 102%).

SUMMARY AND CONCLUSION

A new method was established for simultaneous estimation of Kanamycin and Cephalexin by RP-HPLC method. The chromatographic conditions were successfully developed for the separation of Kanamycin and Cephalexin by using SYMMETRY C18 column $(4.6\times150\text{mm})5\mu$, flow rate was 1ml/min, mobile phase ratio was (70:30 v/v) methanol: phosphate buffer $(KH_2PO_4\text{and }K_2HPO_4)$ pH 3

(pH was adjusted with orthophosphoric acid), detection The wavelength was 258nm. instrument used was WATERS HPLC Auto Sampler, Separation module 2695, photo diode array detector 996, Empower-software version-2. The retention times were found to be 2.403 mins and 3.907 mins. The % purity of Kanamycin and Cephalexin was found to be 100.27% and 99.87% respectively. The system suitability parameters for Kanamycin and Cephalexin such as theoretical plates and tailing factor were found to be 2294, 1.27 and 4891 and 1.03, the resolution was found to be 8.67. The analytical method was validated according to ICH guidelines (ICH, Q2 (R1)). The linearity study n Kanamycin and Cephalexin was found in concentration range of 50µg-250µg and 5µg-50µg.

REFERENCES

- 1. Herman G, Bergman A, Liu F, Stevens C, Wang A, Zeng W, Chen L, Snyder K, Hilliard D, Tanen M, Tanaka W, Meehan A, Lasseter K, Dilzer S, Blum R, Wagner J. Pharmacokinetics and pharmacodynamic effects of the oral DPP-4 inhibitor sitagliptin in middle-aged obese subjects. J.Clin.Pharmacol.2006; 46 (8): 876–886.
- 2. *K. G. Patel et al*, High-Performance Liquid Chromatographic Determination of Cefalexin Monohydrate and Kanamycin Monosulfate with Precolumn Derivatization. Journal of Liquid Chromatography & Related Technologies Volume 38, Issue 6, 2015.
- 3. *Patel satish.A et al*, Spectrophotometric Determination Of Cephalexin Using Ninhydrin Reagent In Tablet Dosage Form. Patel satish A et al.IRJP 2011,2(9),123-126.
- 4. *A. Vyasa et al* ,Development and Validation of Spectrophotometric Method for Estimation of Cephalexin in Bulk and Tablet Dosage Forms. Oriental Journal of Chemistry 2011, Vol. 27, No. (1): Pg. 359-362.
- 5. *Nandan V. Shinde et al*, Development and Validation of HPTLC Method for Estimation of Cephalexin in Human Plasma. International Journal of Research in Pharmaceutical and Biomedical Sciences ISSN: 2229-3701.
- 6. Sagar Suman PANDA et al., Determination of Cephalexin Monohydrate in Pharmaceutical Dosage Form by Stability-Indicating RP-UFLC and UV Spectroscopic Methods. Scientia pharmaceutica 2013; 81: 1029–1041
- 7. *J Vijay et al.*, Formulation And Evaluation Of Cephalexin Extended-Release Matrix Tablets Using Hydroxy Propyl Methyl Cellulose As Rate-Controlling Polymer. J Young Pharm. 2012 Jan-Mar; 4(1): 3–12.

- 8. *Kalpana G. Patel et al.*, Derivative and ratio derivative spectroscopic methods for simultaneous determination of Cefalexin monohydrate and Kanamycin monosulphate using phenylisocyanate as derivatizing reagent. Kalpana G. Patel et al. / Journal of Pharmacy Research 2014,8(5),689-695.
- 9. *Kaale E et al.*, Development And Validation Of A Simple Capillary Zone Electrophoresis Method For The Analysis Of Kanamycin Sulfate With UV Detection After Pre-Capillary Derivatization. J Chromatogr A. 2001 Jul 27;924(1-2):451-8.
- 10. Tripathi K.D. Essential of Medical Pharmacology, 5th Edn, Jaypee Brothers Medical publisher New Delhi. pp: 515-516.
- 11. Inoue S, Bian K, Okamura T, Okunishi H, Toda N (July 1989). "Mechanisms of action of eperisone on isolated dog saphenous arteries and veins" *Jpn. J. Pharmacol.* 50 (3): pg 271–282.